


G 500i

Electrolyte Analyzer



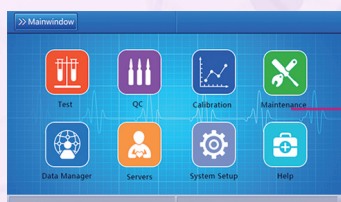

Operation Interface

Thermal Printer




Flow Path

Reagent Pack
and TCO₂ Module



Probe

Electrodes Group

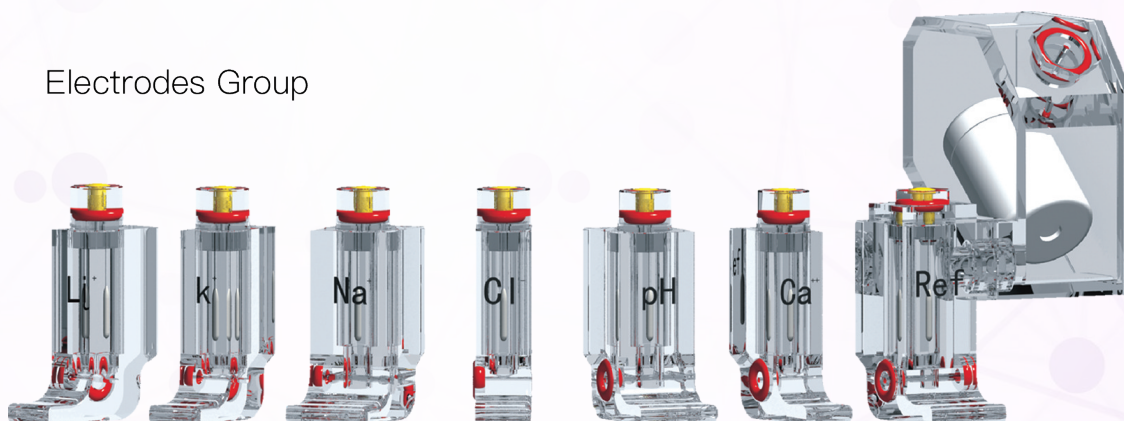

Accuracy and Reliability

- Focus on electrolytes development and innovation for more than 20 years
- Maintenance free electrodes with long shelf life
- Imported TCO₂ module to ensure accuracy and reliability
- Build-in integrated reagent pack to reduce biological contamination
- Auto calibration and slope-offset adjustment to ensure accurate result
- Excellent measuring range and detection precision

Simplicity and Efficiency

- User friendly interface and online fault alarm and troubleshooting
- One-button maintenance procedure to improve efficiency
- Sample volume lower to 80µl, Measuring cycle less than 25 seconds
- Results storage >10,000
- Standard network configuration and online upgradable function
- Barcode management and reagent volume alert

Practical and Economical

- 7 inch colorful touch screen
- Build-in thermal printer
- Upgradable autoloader
- Air and liquid mixture flush to eliminate carryover
- Smart sensor to diagnostic and adjust hardware
- Competitive cost per test

Reagent Pack



Barcode Scanner



Model:

G 500i : K⁺、Na⁺、Cl⁻、Ca²⁺、pH、

Sample Type: Serum, Whole Blood, Plasma, Urine, CSF

Test Range and Accuracy

Item	Range (mmol/L)	SD (mmol/L)
K ⁺	0.5~15.0 mmol/L	≤1.0%
Na ⁺	20.0~200.0 mmol/L	≤1.0%
Cl ⁻	20.0~200.0 mmol/L	≤1.0%
iCa ²⁺	0.1~5.0 mmol/L	≤1.5%
pH	6.0~9.0	≤1.0%
Li ⁺	0.0~3.0 mmol/L	≤1.5%
TCO ₂	6.0~50.0 mmol/L	≤3.0%

Working Environment

Ambient Temperature: 10°C~40°C	Relative Humidity: ≤80%
Atmospheric Pressure: 86kPa~106kPa	Power Supply: 100V~240V~ 50Hz/60Hz
Power: 300VA	Properly grounded, away from electromagnetic interference, avoid direct sunlight
Dimension: 340x283x455mm	Net Weight: about 12kg

Reagent Kit

Item	Model	Packing Specification	Type
ISE Reagent Pack	G 500i-A8 II	One pack including Cal A 1000ml, Cal B 250ml and TCO ₂ calibration 50ml	Liquid
ISE Reagent Pack	G 500i-A8 III	One pack including Cal A 1000ml and Cal B 250ml	Liquid
TCO ₂ Reaction Solution	G 500i-A7 I	One bottle of 100ml	Liquid

All rights reserved. Specifications subject to change without prior notice

Version 2.0



Property No. 14, Second Floor, Industrial Area, Tilak Nagar, New Delhi (India) 110018
 +91-9643001230, 9643001224, 9643001225
 hello@goodwilldiagnostics.in
 www.goodwilldiagnostics.in



ISO9001 & ISO13485